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Description

This test kit allows you to quickly and reliably determine the iron content of your water. The kit is immediately ready for use. After the measurement has been carried out, you can dispose of the sample easily with the domestic waste water.

Function

The analysis is carried out according to the colorimetric method. By adding the reagents, a colour change is made to the parameter to be analysed. You can determine the measurement result by comparing it with a colour chart consisting of 6 colour fields.

Features

- easy to use
- fast and reliable results
- use of high quality chemicals

Remark

Our test kits are produced entirely in Germany. Regular quality control in our in-house laboratory ensures the highest quality of our products. Related documents can be found on our website.

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
Instructions

Please note that the displayed colour may differ from the original colour chart on the screen or a printed version. Always use the chart supplied with the test kit. If you are missing this, we will be pleased to send you a new chart.

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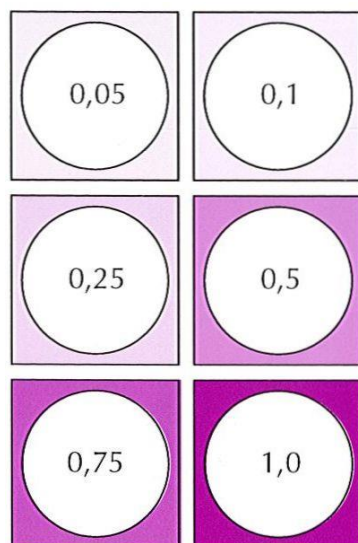Instructions for use:

1. Fill the tube up to the 10 ml-mark with water sample.
2. Add 4 drops of reagent A to the water sample, mix and allow to stand for 3 min.
3. Add 4 drops of reagent B to the water sample, mix and allow to stand for 5 min.
4. Put the tube on the color field and match the color by looking from the top through the tube.